


# Adding a tilt feature to the seat

## TILDA™

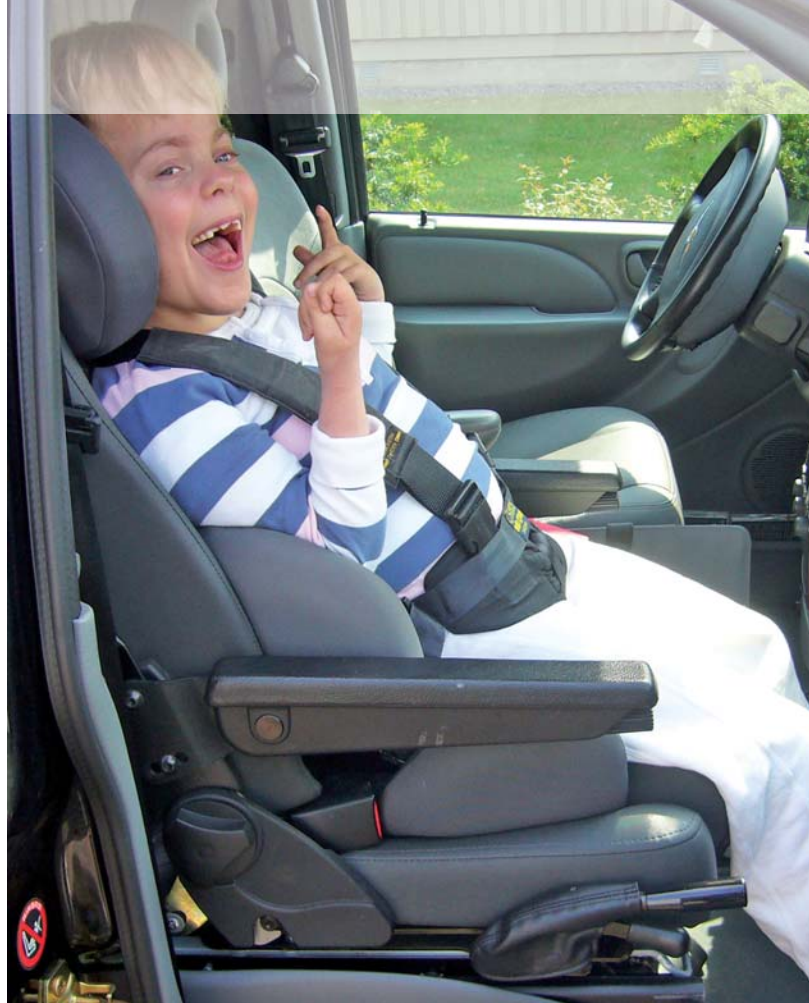

# INVENTIONS FOR QUALITY OF LIFE

# TILDA™

## - tilts the seat in 3 positions

TILDA provides the possibility to tilt the entire seat up to 15°, in three steps.

When seated in tilt position 1, the car ride becomes more relaxing and the adjustability also provides relief when suffering from pain or spastic conditions.

In combination with CARONY, TILDA will provide more head space as the seat slides over from the Carony wheel unit to the swivel base inside the car.

Maximum user weight: 150 kg.  
Height, under the seat: 24 mm Width: 520 mm



*Tilda in position 3 and in combination with the TURNOUT swivel base, with the Carony rails mounted on the sides. The footrest is an option.*

TILDA is specially designed to be combined with:

- TURNOUT swivel base that swivels the seat out over the sill.
- TURNOUT + CARONY which enables the seat to be placed on a wheel unit and become a wheelchair.
- GLIDE RAILS (Scheelman), providing length adjustability for the seat.



*Grab hold of both levers to release the position of the seat, then raise the seat to an upright position.*



## Crash tested

When tilted in position 1 the ride will become more comfortable and the sitting more relaxed. TILDA is crash tested in pos. 0 (upright position) and in pos. 1

TILDA is in compliance with following regulations and directives:

74/60/EEC - 2000/4/EG, 74/408/EEC - 2005/39/EEC  
93/42/EEC, 95/28/EG, ECE R-17, ECE R- 21,  
SS-EN 12182, ISO 3795  
FMVSS 201, 207, 208 and 302